




**Schlein**  
**Ultra™**  
Shoulder Positioner





*OSI's Third Generation Shoulder Positioner...*

# Schlein Ultra<sup>™</sup> Shoulder Positioner

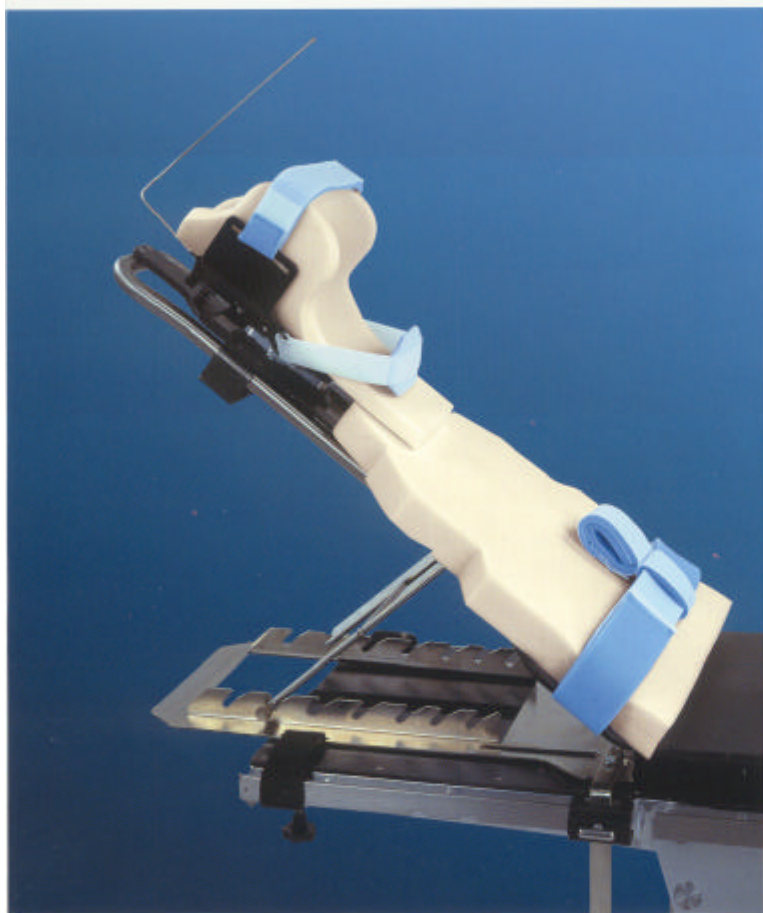


The Schlein Ultra Shoulder Positioner provides a safe, convenient and stable method of maintaining patients in the beach chair position for arthroscopic and open shoulder surgery. The setup is minimal and concerns associated with lateral decubitus positioning are eliminated. The Schlein Ultra fits all operating room tables and offers an effective solution for shoulder surgery positioning requirements.



Headrest wings retract posteriorly creating a totally obstruction-free space around the patient's head for the anesthetist.

With the wings in place and head strap positioned, the headrest becomes a secure head support.

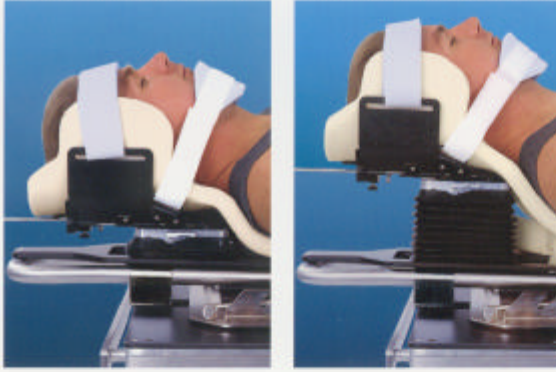

The cervical flexion angle for a kyphotic patient is accommodated with the adjustable headrest.

The Schlein Ultra is designed to suit the high demands of the OR Suite.

- The center mounted support rod assures unobstructed radiolucency at the surgical site.
- The phenolic backrest has a low contour for optimal shoulder access both anteriorly and posteriorly.
- All release levers and adjustment knobs are readily located for convenient access and easier positioning of the patient.
- The headrest adjusts superiorly and inferiorly to accommodate the patient's height.
- Single use Patient Care Kits may be customized to accommodate each patient. On the operative side, tear along the perforations to increase shoulder access. Leave pad in place to support the non-operative shoulder.



The headrest tilts laterally up to 30° in either direction.

## ORDERING INFORMATION



**5358  
Schlein Ultra  
Shoulder Positioner**  
Comes with three Patient  
Care Kits (OSI Model No. 5348).

**5345  
Patient Care Kit**  
Three per case.  
(Others available.)

**5394  
Aluminum Side  
Rail Socket**  
(Two required.)




**SUSP  
Schlein Ultra  
Complete**

**5338-24  
Schlein Shoulder  
Positioner Cart**

This specially designed  
storage cart transports the  
positioner, accessories and  
reserve Patient Care Kits.

*(The Schlein Ultra fits  
in the same storage  
cart for the SSPIL.)*

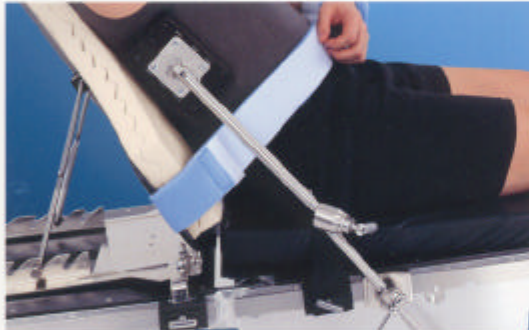

### **5338-100 Lateral Support Pad Assembly**

This specialty pad assembly provides additional support and prevents the patient from shifting laterally during the rigors of surgery. (Clark Socket included.)



### **5338-70 Translating Arm Board for Schlein**

Arm board may be positioned to provide a platform for either arm.

### **5393 Clark Socket**

## SPECIFICATIONS

Headrest Translation	5"
Anterior Range	8.5"
Longitudinal Range	±30°
Lateral Tilt	400 lbs.
Patient Weight Limit	29 lbs.
Schlein Ultra Weight	

*Attaches to any OR table.*

## WARRANTY

One year parts and labor against defects in materials or manufacturing.

The manufacturer thanks Allan P. Schlein, M.D., for his assistance in the development of the OSI Schlein Shoulder Positioner.

U.S. Patent No. 5,515,897. Other patents pending.